

User Manual

Shower stool

Neapel	
REF	368.00.40
	130 kg



Figure similar



Read this User Manual **before using the product for the first time** and keep it at hand for the user as well as for future use and provide it to the purchaser in the case of subsequent sales.



Table of contents

1. Product description	3
1.1. Introduction	3
2. Features	3
2.1. Intended purpose	3
2.2. Use	3
2.3. Storage & transport.....	4
2.4. Equipment & functions	4
2.5. Packaging	4
2.6. Disposal	4
3. Symbols & nameplates	4
4. Safety instructions	4
5. Incident reporting.....	5
6. Warranty.....	5
6.1. Complaints.....	5
7. Reuse, lifetime	5
8. Scope of delivery.....	5
9. Setting up for use.....	6
9.1. Assembling the Tube Sets	6
9.2. Adjusting the seat height	6
10. Product information	6
11. Operation of the product	6
11.1. Sitting down	7
11.2. Getting up.....	7
12. Cleaning & disinfection.....	7
12.1. General	7
13. Maintenance.....	7
14. Spare parts & accessories.....	8
15. Drawing & dimensions.....	9

1. Product description

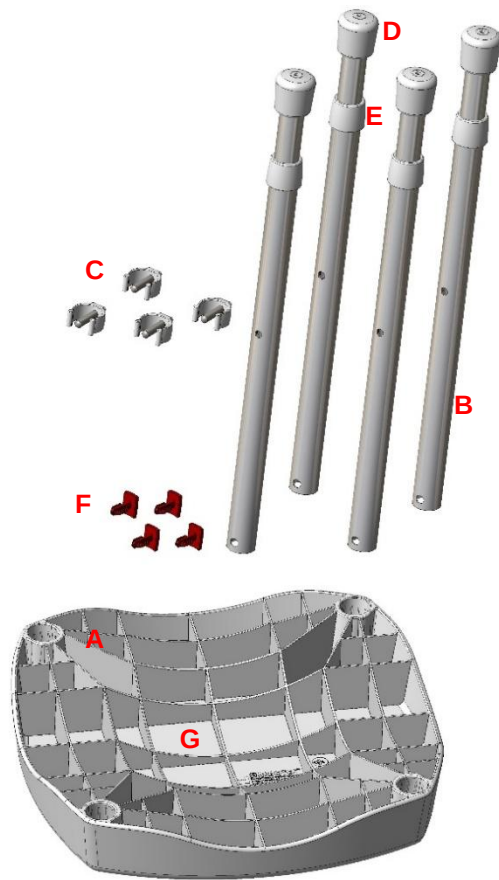


Figure shows:

A	Seat plate
B	Tube set
C	Adjustment clip
D	Foot cap
E	Guide bush
F	Pine-tree clip
G	Nameplate

1.1. Introduction

- The product described here is a medical device of risk class I according to (EU) 2017/745 MDR.
- In this user manual, you will find information on the safe and proper initial operation, use, maintenance and care of this product.
- Read this manual carefully before using this product. Follow the safety instructions and act accordingly.

2. Features

2.1. Intended purpose

- The product is intended for the care of patients with limited ability to stand. It allows you to shower while sitting down. The patient must have sufficient sitting ability and upper body musculature. If necessary, a caregiver must be present during use. The product is not intended for prolonged sitting.

Specification:

- Maximum user weight: 130 kg
- Stool
- Indoor use

Indication

- If it is impossible for patients to safely stand in the shower, due to damage to the function of the lower extremities with limited or eliminated mobility and muscle strength and/or damage to the vestibular functions or tendency to fall due to other causes, while at the same time maintaining a sufficiently stable thorax posture for sitting and sufficient residual functions of the upper extremities.
- Purpose of care: Enabling independent or more independent showering

Example diagnoses:

- The spectrum of possible diagnoses is extensive. It ranges from general weakness due to old age to arthrosis of the lower extremities.

Contraindication

- Functional impairments of the lower extremities, trunk posture, or vestibular functions that preclude the safe use of the shower stool. Balance disorders and general sensory impairments that could compromise safety must be assessed by qualified professionals.
- The product must only be used on intact skin.

2.2. Use

- Temperature range: 5 °C to +35 °C Direct radiant heat on the product must be avoided.
- Relative humidity: 40%-75%. Outside of these conditions, the product must only be used for a short time.
- The product is not suitable for use in swimming pools or saunas, as surface discolouration and limitations of mechanical strength may occur.
- Disassembly for transport or storage is not possible.

- There must be a flat surface of sufficient size. See 15.

2.3. Storage & transport

- Temperature range: 5 °C to +25 °C; direct radiant heat on the product must be avoided.
- Store the product in dry rooms with a humidity of up to 65%.

2.4. Equipment & functions

- Adjustable height with 6 positions

2.5. Packaging

- We recommend you keep the packaging for later transport.
- Dimensions: 363 mm x 367 mm x 55 mm

2.6. Disposal

- Disposal and recycling of used products and packaging must be carried out in accordance with currently valid regulations. Please contact a disposal company for further information.
- The Infection Protection Act must be observed.
- Please observe any labels and information on packaging material and act accordingly.

3. Symbols & nameplates

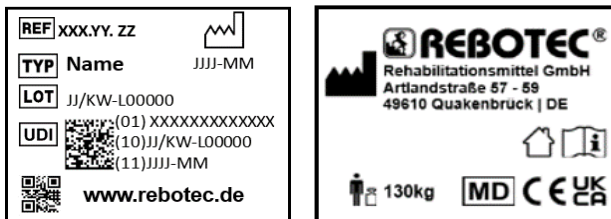


Figure similar.



The nameplate must not be removed! For ordering spare parts and processing complaints, the device name, serial number and year of manufacture will be required.

Symbol	Bedeutung	Symbol	Bedeutung
	Product number		Item description / Proper name
	Batch number		Serial number
	Unique identifier of a medical device		Medical device
	Important information!		Attention! Observe safety instructions.
	Language		Maximum user weight in kg
	For indoor use only.		Suitable for indoor and outdoor use
	Manufacturer		date of manufacture
	Read the operating instructions!		Follow the instructions for use.
	Product complies with (EU) 2017/745 MDR		Product complies with UK MDR 2002
	Protect from sunlight		Store in a dry place
	Temperature limitation		Humidity, limitation
	This page at the top		Fragile, handle with care
	Care instructions		Maintenance notice
	Do not dispose of in household waste.		Stack limit n = number

4. Safety instructions



Please observe the following safety instructions:

- Before using the product for the first time, the User Manual must be read completely in order to avoid hazards or damage caused by maloperation.
- Before every commissioning, check whether there is any visible damage to the product.
- In case of malfunction or damage of any kind, the product must not be used.
- Do not leave children unattended with the product.

Caution: Danger of falling!

- Be sure to check the following before use:
 - Stability of the product.
 - Seat plate for damage.
 - Adjustment clips of the height adjustment for tight fit.
- Do not carry the product together with the patient.
- Only use the product indoors and on a flat surface.
- Only use the product if completely assembled.
- Do not slump carelessly on the product.
- The user must sit completely on the seat.
- Avoid excessively leaning over the product.
- Avoid wobbling or rocking the product.
- During use, all legs of the product must be in contact with the floor.
- Adjust all four sets of tubes to the same height.
- Only use the height adjustment in the holes provided for it.
- Attention! Never use the product without the adjustment clips in place.

Caution: Danger of breakage!

- The product is only suitable as a seat for persons. Observe the permitted user weight (2.1).
- Do not expose the product to strong temperature fluctuations (2.2 & 2.3).

Caution: Risk of injury!

- There is a risk of injury with surface temperatures > 40 °C. Therefore, do not expose the product to high temperatures (e.g. sunlight, radiators, hot water). Allow the product to cool down before use.



If the product is used by patients with insensitive skin (no temperature perception or skin damage), a heat check (e.g. touching with the back of the hand) must be carried out by the nursing staff.

5. Incident reporting

- Serious incidents related to the product must be reported immediately to the manufacturer and the competent authority.
 - Use the reporting functions and provided forms.
- Manufacturer: pms@rebotec.de

6. Warranty

- The manufacturer grants a warranty of 12 months for this product. The prerequisites are the terms and conditions (www.rebotec.de/agbs), as well as the intended use. The statutory warranty provisions shall apply.
- Unauthorised modifications to this product will void the product's conformity and the warranty.

6.1. Complaints

- Please contact us before returning the product.
- To reduce transport damage, use the original packaging if possible.
- The Infection Protection Act must be complied with.
 - Please make sure that the product does not pose a risk of infection when returned.
 - Enclose the information sheet on safety with the product.
 - Costs may be incurred for returns that are not marked as harmless.

7. Reuse, lifetime

- The product is suitable for reuse. How often the product can be reused depends on the frequency and nature of use.
- Before reuse, the product must be reprocessed hygienically and the technical and functional safety and reliability of the product must be checked and restored if necessary. See 12.
 - For more detailed information, please see the reprocessing instructions (02.12.159) and the maintenance instructions (02.11.163).
 - The documents can be downloaded from the following link.
 - www.rebotec.de/downloads
 - A maintenance record is included in the service manual to document service and repair.
- If the product is used as intended, its lifetime is 3 years. The effective lifetime may vary depending on the frequency and intensity of use and the general condition. A longer service life is possible with a lower intensity of use.

8. Scope of delivery

- 1x Shower stool Neapel
 - 1x seat plate
 - 4x tube set


- o 5x Pine-tree clip
- 1x User Manual (02.08.212)
- The reprocessing and maintenance instructions can be ordered or viewed and downloaded at www.rebotec.de/downloads.

9. Setting up for use

- The product has been carefully inspected in the factory for freedom from defects and completeness. Please check the product for any possible damage which may have occurred during transport and for completeness of the scope of delivery upon receipt.
- Observe the separate instructions for assembly (02.13.211).

9.1. Assembling the Tube Sets


- Place the seat plate with the underside facing up on a clean and level surface.
 - o Use a soft surface if necessary to avoid scratches.

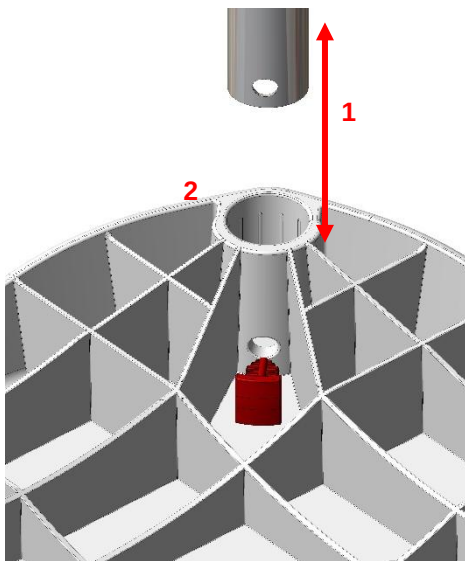

- Slide the tube set (1) into the socket (2) until it is fully inserted.
- After inserting the tube assembly, insert the Pine-tree clip into the tube until it is fully seated and secure.
- The Pine-tree clip secures the tube in the mount and prevents it from coming loose accidentally.
- Once all tube assemblies are installed, turn the item over.
 - o No tube assembly should fall out.
 - o When correctly assembled, the item must not wobble on a flat surface.

9.2. Adjusting the seat height

- Place the seat plate on a clean and flat surface with the underside facing upwards.

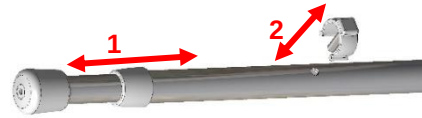

- Remove the clip (2).
- Pull out the inner tube (1) to the desired position.
 - o The spacing between the individual positions is 25 mm.
- Re-insert the clip (2).
 - o The adjustment clip must always engage in one of the holes of the inner tube.
 - o Never use the product without the adjustment clips in place.
 - o All tube sets must be set to the same height.
- All four legs must be set to the same height.



The item is not designed to be disassembled. The legs are fixed to the seat panel in a way that does not allow for non-destructive removal.

The article is not suitable as a travel stool!

10. Product information

- We reserve the right to make technical developments and improvements with respect to the design shown in this user manual.
- See 15

Field	Value
User weight	max. 130 kg
*Product weight	1.486 kg
Total weight	1.692 kg
Packaging unit 12 pieces	21,597kg
Skin sensitivity	DIN EN 10993-1/-5
Flammability**	DIN EN 1021-1/-2

**Without accessories.*

***Flame retardant: For products made of PUR material.*

11. Operation of the product



Attention: Danger of falling!

- Observe the safety instruction (4).
- Before use, it must be ensured that the product is in proper condition.

- A stable position with no risk of slipping or tipping over must be ensured.

11.1. Sitting down

- The user must position themselves as close as possible to the product. If necessary, he/she can orientate or support himself/herself with his/her hands on the edge of the seat plate.
- If this is not possible due to the patient's physical limitations, the user must seek the assistance of a competent person. This can be done with the help of an assistive device.

11.2. Getting up

- To get up, the feet must be in contact with the floor.
- Then, slightly shift your weight forward and slowly rise from the seat.
- Do not lean your body too far forward to avoid the risk of falling. If necessary, the user must seek help to get up from the stool.

12. Cleaning & disinfection

12.1. General

- The item is suitable for manual cleaning and disinfection.
- Use the validated agents in accordance with our reprocessing instructions.
- Observe the reprocessing instructions (Doc. no.: 01.12.159).
- You can order the reprocessing instructions or view and download them at
 - www.rebotec.de/downloads

13. Maintenance

- Wear and tear that is not recognised in time and/or improper handling, as well as irregular or no maintenance may result in accidents.
- The product is maintenance-free when used as intended.
- Regular checks are required:
 - Check the product and accessories for visible damage (cracks, breaks, loose parts)
 - Check stability before use
 - Function of the height adjustment mechanism and clips
 - Locking of the tubes in the seat plate

- If you have any questions concerning the maintenance of the product, please contact an authorised dealer or a medical and healthcare supplies store. In the event of damage or a need for repairs, contact your specialist dealer immediately.


- Do **not** make **any** unauthorised changes to the article. This may affect the safety and function of the item. **In such cases, REBOTEC will not assume any liability.**

- Only use the manufacturer's spare parts for repairs. See 14.

- After maintenance or repair work has been carried out, we recommend you clean the product before using it again.

- Please observe the maintenance instructions (02.11.163).

14. Spare parts & accessories

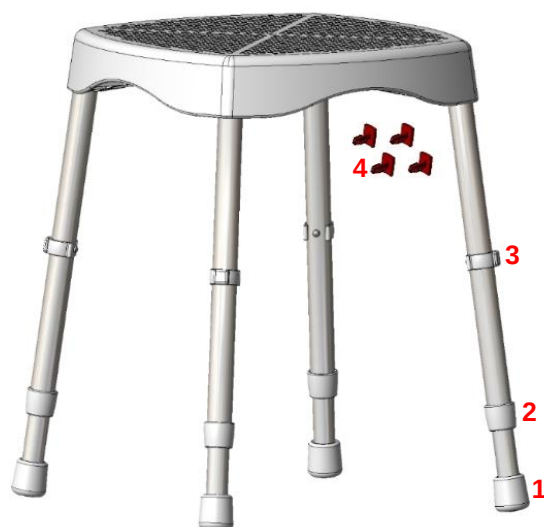


No.	Prod. no.	Designation
1	468.00.19	Foot cap
2	202.90.10	Guide bush
3	201.90.10	Adjustment clip
4	468.44.50	Pine-tree clip
-	02.08.212	User Manual Neapel
-	02.12.159	Reprocessing Instructions
-	02.11.163	Maintenance Instructions



Use only REBOTEC spare parts and accessories to ensure the proper and safe functioning of the product.

15. Drawing & dimensions

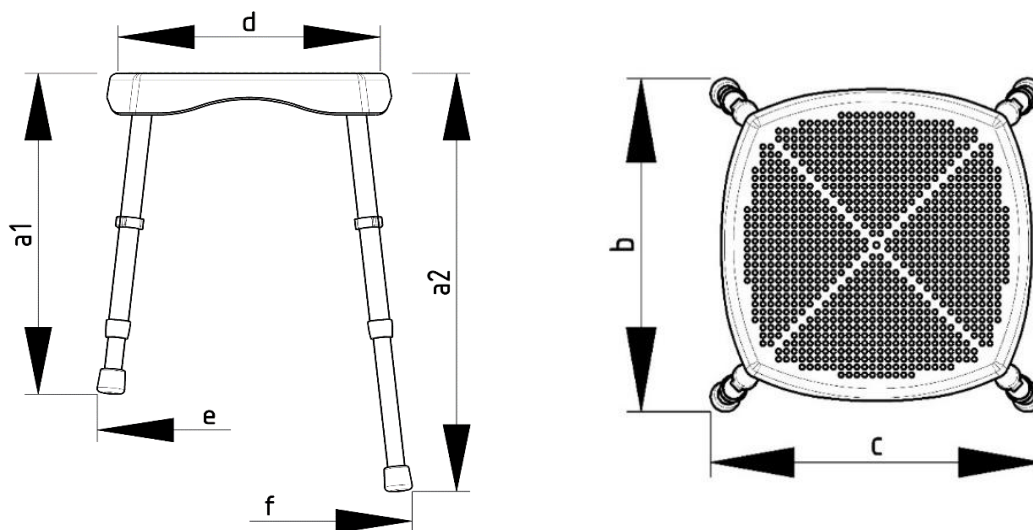

No.	Designation	Dimensions in mm
a1	Height min.	410
a2	Height max.	525
b	Depth	see e & f
c	Width	see e & f
d	Seat width	290
	Seat depth	290
e	Floor space min.	340 x 340
f	Floor space max.	370 x 370



Dimensions and weights may vary depending on the equipment. There might be minor deviations due to production tolerances. Subject to change without notice.





REBOTEC® Rehabilitationsmittel GmbH
D-49610 Quakenbrück, Artlandstr. 57- 59
Telephone: +49 (0) 5431/9416-
Fax: +49 (0) 5431/9416- 66
Web: www.rebotec.de
Email: info@rebotec.de






EN




REBOTEC®

1.

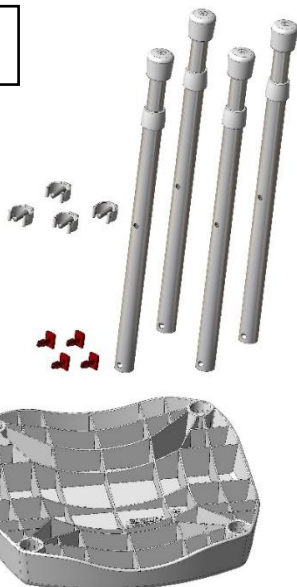


2.

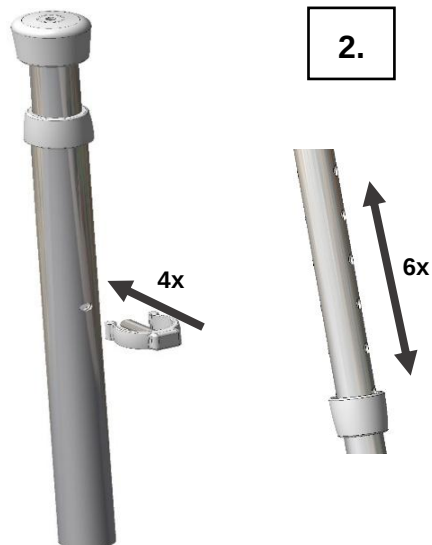

3.




4.

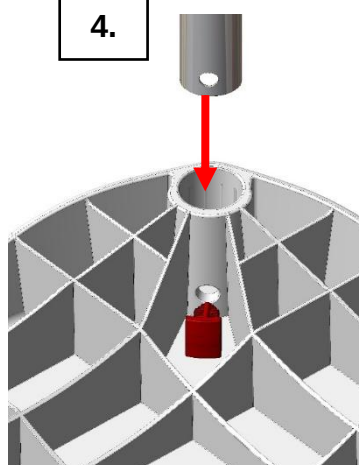


5.



6.

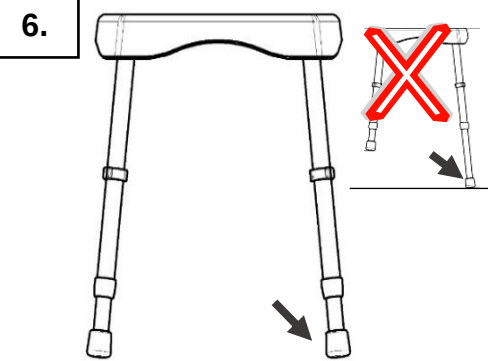

7.

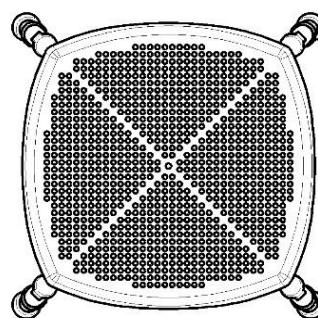

Read User Manual **before using the product for the first time** and keep it at hand for the user as well as for future use and provide it to the purchaser in case of subsequent sales.